

EDGE 3



The Edge 3 is the culmination of 8 years of product improvement and innovation. Upgraded SRS (Smooth Ride Suspension) provides a smooth, comfortable ride and advanced stability. Optional iLevel® is available at 4.5 mph with up to 12" of power adjustable seat height. Standard features include LED fender lights, a USB charger, and drive wheel colour accents. The Edge 3 is available in a total of 14 colour options, including seven brand-new, bold colours.

UNIQUE TECHNOLOGY

to the industry.

12" iLevel™ | SRS Suspension | TB3® Seating System
Q-Logic 3 with Bluetooth | TB3 Flip-up armrests



QUANTUM

Dedicated to revolutionizing the power chair experience.



CRASH
TESTED



@QUANTUMREHABUK





SRS suspension

A smoother ride experience. SRS technology is automotive-grade, dampened suspension making for a more comfortable ride.

55 Ah Batteries

More miles per charge! 55 Ah batteries allow for increased range, for the most active user. Easy access through front shroud.

New Colours

8 new bold colours to select from including our new carbon fibre shrouds.

Turning Radius

Access smaller spaces easier with a tight turning radius of 20.5" (52.07 cm)

Motor Package

2-Pole advanced responsive motors. 4 Mph or 6 Mph motor option
Accu-trac motor option delivering enhanced traction and tracking

Measurements

14" Drive wheel, 6" low impact solid castor wheels front and rear.
Length: 32.23" (Without legrests) | Width: 24.4" |
Ground clearance: 2.7"

Seating Solutions



TB-Flex (Coming Soon)

- 15" x 12" to 21" x 20" seat size options
- Visco padded or tension adjustable backrest options
- Cost effective power options
- Durable and easily adjustable
- Easy centre of gravity adjustment
- Accepts a wide variety of after market seats and backs



TRU-Balance® 3 Seating

- 12" x 12" to 22" x 22" seat size options
- Comprehensive range of power options
- Available with 12" iLevel Seat elevation
- Accepts a wide variety of after market seats and backs
- Endomorph and mesomorph seating, allows for independent adjustment of seat width and backrest width
- Premium accessories available

Power Options



TRU-Balance® 3 Power Tilt

- 0-50° of tilt
- Lower standard seat-to-floor height
- Easy adjustment and durability
- Maintains centre of gravity shift



TRU-Balance® 3 Power Recline

- 85° - 175° of recline
- Adjustable shear reduction
- Compatible with after-market backrest interface

Powered Elevating leg rests

- Articulating foot platform
- Articulating swing-away leg rests
- Can be operated in conjunction with TB3 recline




12" iLevel™ Power Adjustable Seat Lift

- Lower standard seat-to-floor height
- iLevel™ allows 12" of seat elevation while at a speed of up to 4.5 mph.
- iLevel is retrofittable on the Edge 3



Additional Features

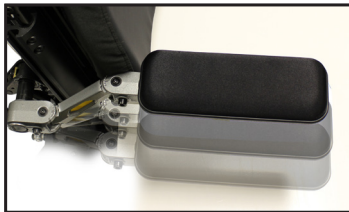

TRU-Balance® 3 Cantilever Flip-up armrests

- Quick height adjustment
- Can be internally or externally rotated for optimal positioning
- Multiple armrest pad options



TRU-Balance® 3 Centre mount foot platform

- Fixed or powered articulating
- Modular footplate extension kits
- Multi-adjustable calf pads
- Paediatric & adult footplate options



Swing-away legrests

- 60° - 90° standard swing-away
- Paediatric & adult footplate options
- Powered Articulating options
- Manual elevating



Attendant/Dual control options

- Full proportional joystick
- Quick and easy to use control interface
- Control the chair at the attendants speed




Body positioning

- A full range of multi adjustable swing-away laterals, offering advanced positioning of the trunk
- Designed to easily allow head to toe positioning using any of the Stealth positioning accessories



After-market backrest interface

- Allows for an unprecedented level of compatibility between TRU-Balance 3 back rest and most after-market backrests
- Retains full functionality, including shear reduction

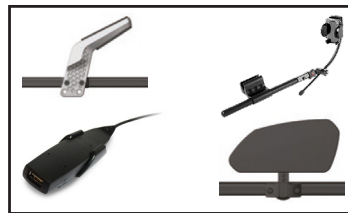

Premium accessories

- Transfer bars
- Mobile phone holder
- Portable USB charger
- Clothes guard
- Additional accessories available



Full range of control devices available

- Switches
- Sensors
- Proportionals
- i-Drive™



OPT: Extended Traction Tyres⁴

Edge 3 | Colour Options & Tech Spec

Colour Options



Cherry Bomb



Back in Black



Sugar Plumb



Green Machine



Grey Street



Evolution Orange



Night Train



Blue Jean



White Satin



Carbon Fibre



Military Green



Raspberry Berry



Lemon Crush



Midnight Blue

Technical Specifications

Drive Wheels (cm) 35.56 cm pneumatic or solid	Caster Wheels (cm) Front and Rear: 15.24 cm	Max Speed (mph) 4 mph or 6 mph with lights ¹	Ground Clearance (cm) Battery Tray: 6.9 cm	Turning Radius (cm) 52.1 cm (20.5")	Overall Length (cm) Without legrests 81.86 cm (32.23")	Overall Width (cm) 61.97 cm (24.4")	Seating Options TB Flex and TB3 Seating	Batteries 2 x 55 AH Deep Cycle	Available Electronics NE+ / QLogic 3	Weight Capacity 21 Stone 6 lbs (136 k.g)
Seat to ground height (") 16.5"-19.5"	Base Weight (kg) 62.86 kg	Battery Weight Each (kg) 17.24 kg	Warranty <ol style="list-style-type: none"> Lifetime limited warranty on frame. 1-year limited warranty on electronics. 18-month warranty on drive control system. 2-year limited warranty on drive train. 			Notes <ol style="list-style-type: none"> Speed varies with user weight, terrain type, battery charge, battery condition and tyre pressure. Due to manufacturing tolerances and continual product improvement, this specification can be subject to variance of + or - 3%. Dimensions listed are for power base only. Overall measurements will vary based on seating and accessory selections. Battery weight may vary based on manufacturer. 				

Get Social with Quantum



FOLLOW THE LIVES OF OUR QUANTUM STARS, AND GET THEIR PERSPECTIVE OF THE EDGE 3

www.quantumrehab.co.uk/quantum-stars

#ileveluk

#myquantumproduct

#socialmobility



@QUANTUMREHABUK

Ensure you have the most up to date information about Quantum Rehab by subscribing to our social media pages

Further information can be obtained through our web site, or prescription forms. All information is subject to change without notification. Please consult Quantum with any queries you may have.